


myQuickmat Prime

new!
be a
flamingo
in a flock
of pigeons

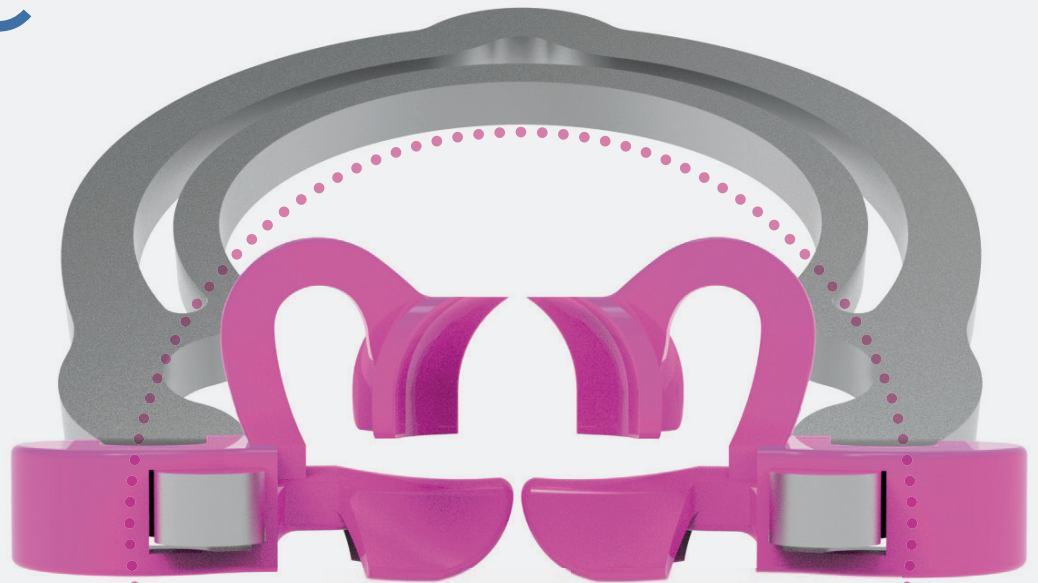

Advantages

simplified posterior system
that reduces chair time

restorative essentials: one
ring, one ring extremity
with integrated wedges and
one anatomical matrix (no
separate wedge required)

self-centring wedge-shaped
extremities: naturally guide
insertion, ensuring effective
bilateral separation of the
adjacent teeth with optimal
matrix adaptation and cervical
seal

ring stability: prevents the ring
from slipping out

matrix contoured from cervical
to occlusal allows accurate
positioning of the contact
point close to the equator

no compromise on efficiency
or performance: the result is
an easier and more natural
restoration of the proximal
anatomy



say
goodbye
to the
wedge



myTines
Prime
flamingo shaped
extremities
with
integrated
wedges

sectional
matrix
system



polydentia
swiss manufacture 



new! myQuickmat Prime



say
goodbye
to the wedge

Testimonial



"The myQuickmat Prime is a revolutionary concept that not only simplifies the protocol for Class II posterior restorations, but also makes the entire process easier and faster. Achieving proper curvature of the proximal contact has always been a challenge for Class II restorations, but with the combination of the new Quickmat Prime and myRing Prime, it has never been easier or more cost effective. In short, a real game changer for my Class II restorations."

Dr. Javier Tapia,
Spain



Images courtesy of Dr. Javier Tapia



for more
information



tines form the matrix
very effectively on
the proximal walls

flamingo shaped

myTines Prime

wedge-shaped reusable
and autoclavable extremities
for myRing Prime



BILATERAL INTEGRATED WEDGES



- self-centring: naturally guide insertion
- effectively separate the adjacent teeth
- form the matrix on the cervical walls
- prevent the ring from slipping out

myRing Prime

stainless steel
sectional
matrix ring

- bent ring so as not to interfere with adjacent teeth



Quickmat Prime

anatomical sectional
matrices, contoured
from cervical to occlusal



- accurate position of the contact point close to the equator
- one size for both premolars and molars
- central upper tab for easy positioning
- rigid but easy to burnish

Quickmat Prime

- natural curved anatomy
- contact point close to the equator

Traditional sectional matrix

- "straight" tooth anatomy
- high contact point



watch
the video

myQuickmat Prime is a sectional matrix system designed to simplify the overall Class II restoration procedure and reduce chair time. The removal of the wedge as part of the restorative procedure makes the system simpler, faster and smarter than traditional sectional matrix systems.

myQuickmat Prime kit REF 6890

- 2 pcs. myRing Prime
- 10 pairs myTines Prime
- 100 pcs. Quickmat Prime matrices



refills



myRing Prime - REF 6891
2 pcs. myRing Prime, 4 pairs myTines Prime



myTines Prime - REF 6892
6 pairs myTines Prime



Quickmat Prime - REF 6893
100 pcs. anatomical sectional matrices